

BRHNY-Series Bio-Burden Reduction Grade Nylon Plus+

BRHNY-Series Bio-Burden Reduction Grade Nylon Plus+ filter cartridges feature Nylon 6,6 membrane with an advanced positively-charged surface modification that is highly efficient in capturing submicronic particulate matter much finer than the stated mechanical rating. This offers a well-proven capability for highly efficient retention of haze, colloids, and color bodies. Specific to its use in medical applications, pyrogenic endotoxins are effectively removed as is well-documented in field use, industry journals, and laboratory data. Superior microbial retention is achieved to deliver a stable and consistent effluent. The BRHNY+ series offers a more cost-effective alternative to hollow-fiber cartridges in many high-purity applications. Cartridges are manufactured in a cleanroom environment and are flushed with 18 megaohm ultra-high purity water to ensure cleanliness, low extractables, and quick rinse-up for service use. Tolerant of repeated hot water sanitization and in-situ steam sterilization cycles for maximum service life. Each element is 100% integrity tested to Global Filter standards to assure consistent and optimal performance.

Construction Materials

Bacterial endotoxin is the pyrogen of greatestconcern in the pharmaceutical and medical device industries. BRHNY+ filter elements have demonstrated capability to remove bacterial endotoxin to below a 0.005 EU/milliliter detection limit at all data points in independent testing.

Microbial Retention Performance

Rating	Challenge Microbe	Log Reduction Value (LVR)		
0.05μ		>10.1		
0.10μ	Brevundimonas diminuta	>9.1		
0.2μ		>9.0		


Typical Applications

- Medical device reprocessing
- Water for Injection (WFI)

Endotoxin removal

Pyrogen removal

Flow Rate vs Pressure Drop

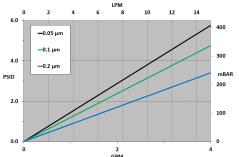

Construction Materials

Membrane	embranePositively-charged Nylon 6,6				
	on polyester substrate				
Support Media	Polypropylene				
End Caps	Polypropylene				
Center Core	Polypropylene				
Outer Support Cages	Polypropylene				
O-Rings/Gaskets	EPDM, Silicone, Viton®				

Dimension (Nominal)

Length	.10 to 40 in (25.4 to 101.6 cm	٦)
Outside Diameter	2.7 in (7.06 cm	۱)

Operating Conditions

Change Out ΔP (recommended).	35 PSID (2.4 bar)
Temperature (max)	176°F (80°C)
Differential Pressure (max)	50 PSID at 68°F
	(3.4 bar at 20°C)

Sanitization/Sterilization

Filtered Hot Water 176°F (80°C) for 30 min
Steaming Sterilization 250°F (121°C) for 30 min
multiple cycles

Chemicals: Cartridges are compatible with most chemical sanitizing agents.

Note: Stainless steel insert option required for all cartridges being hot water sanitized or steam sterilized.

All polypropylene components meet the specifications for biological safety per USP Class VI - 121°C for plastics.

Food Safety Compliance

Materials of construction comply with FDA regulations for food and beverage contact use as detailed in the US Code of Federal Regulations, 21CFR. Materials used to produce filter media and hardware are deemed safe for use in contact with foodstuffs in accordance with EU Directives 1935/2004, and/or 10/2011.

AAMI Standard ST108 Compliance

The BRHNY+ is a key component of water systems required to be compliant to AAMI Standard ST108: Water for the Reprocessing of Medical Devices. The BRHNY+ delivers highly efficient scavenging of microbes and endotoxin when used in a properly-designed system maintained to ST108 recommended practices. This assures continued compliance of the system to meet the highest standards of cleanliness and user confidence.

Ordering Information

•							
BRHNY+	Rating(µ)	Α	Length	С	End Cap Style	O-Rings/Gaskets	Adders
	0.05		10" (25.4cm)		3 = 222 w/Fin	E = EPDM	CS = 316SS Compression Spring
	0.10		20" (50.8 cm)		4 = 222 w/Flat Cap	S = Silicone	I = Stainless Steel Insert
	0.2		30" (76.2 cm)		6 = 226 w/Flat Cap	V = Viton®	
			40" (101.6 cm)		7 = 226 w/Fin		
					16 = 213 Internal O-Ring		
					28 - 222 3-Tahs w/ Fin		

DISCLAIMER: Filtration data presented is representative of performance observed in controlled laboratory testing. It is not given as a warranty, specification or statement of fitness for use. Specific performance can vary widely depending on contaminant type, fluid properties, flow rates and environmental conditions. It is recommended that users conduct thorough qualification testing to assure the product functions as required. For additional technical support, a product Performance Guide is available upon request. SV BRHNY 0524



Email: gfcustomerservice@filtrationgroup.com